
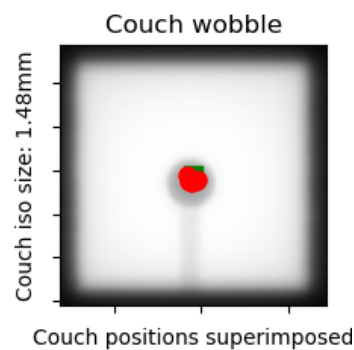
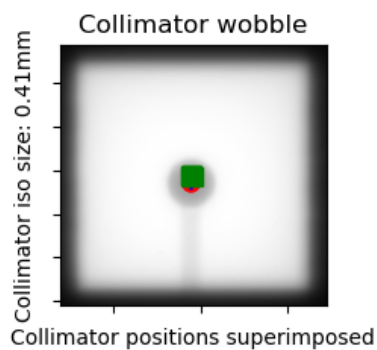
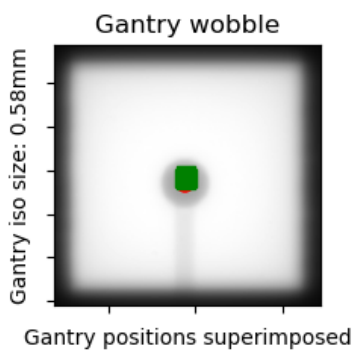
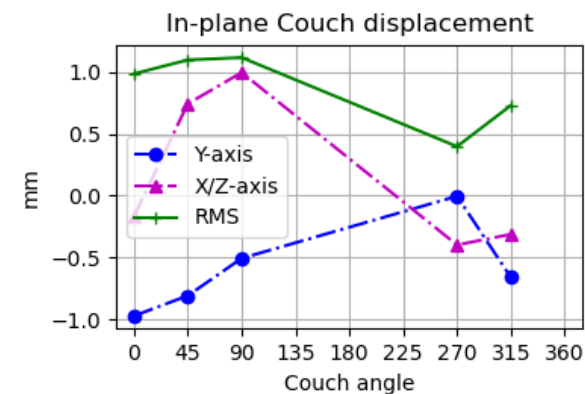
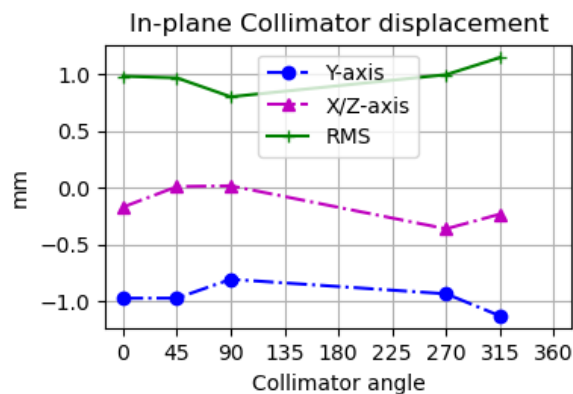
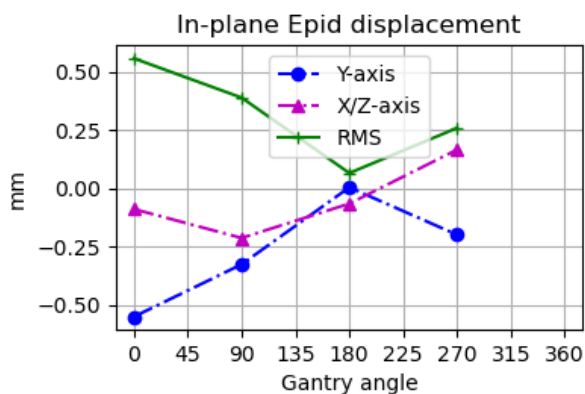



Winston-Lutz Analysis

Winston-Lutz Analysis


=====

Number of images: 12
Maximum 2D CAX->BB distance: 1.15mm
Median 2D CAX->BB distance: 0.89mm
Shift to iso: facing gantry, move BB: RIGHT 0.10mm; IN 0.81mm; UP 0.30mm
Gantry 3D isocenter diameter: 0.58mm (4/12 images considered)
Maximum Gantry RMS deviation (mm): 0.77mm
Maximum EPID RMS deviation (mm): 0.73mm
Gantry+Collimator 3D isocenter diameter: 0.75mm (8/12 images considered)
Collimator 2D isocenter diameter: 0.41mm (5/12 images considered)
Maximum Collimator RMS deviation (mm): 1.15
Couch 2D isocenter diameter: 1.48mm (5/12 images considered)
Maximum Couch RMS deviation (mm): 1.12

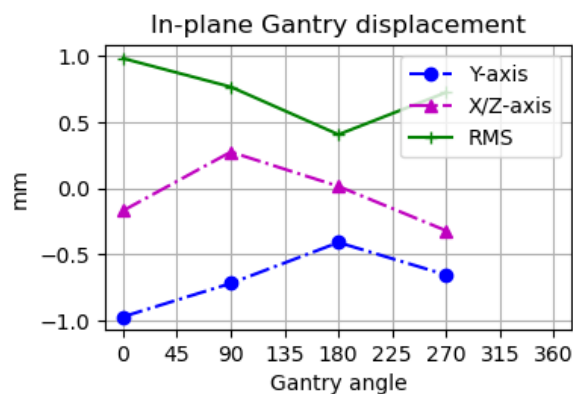

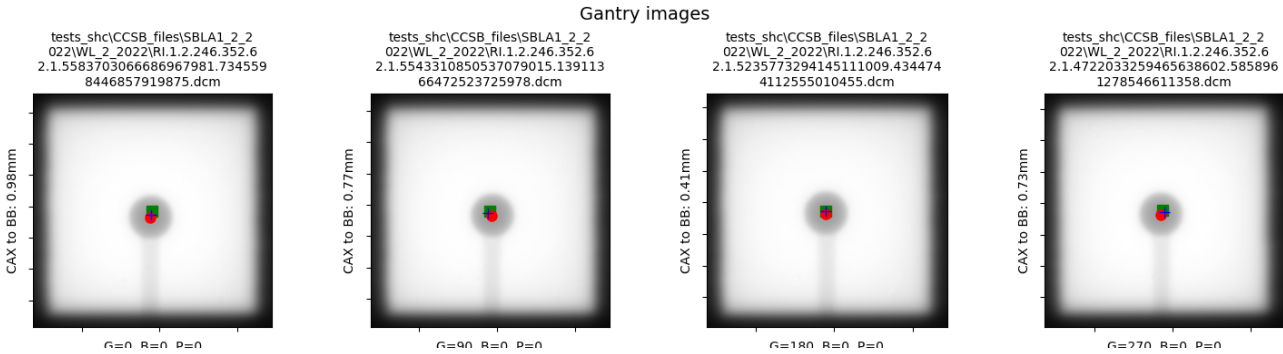

Collimator images

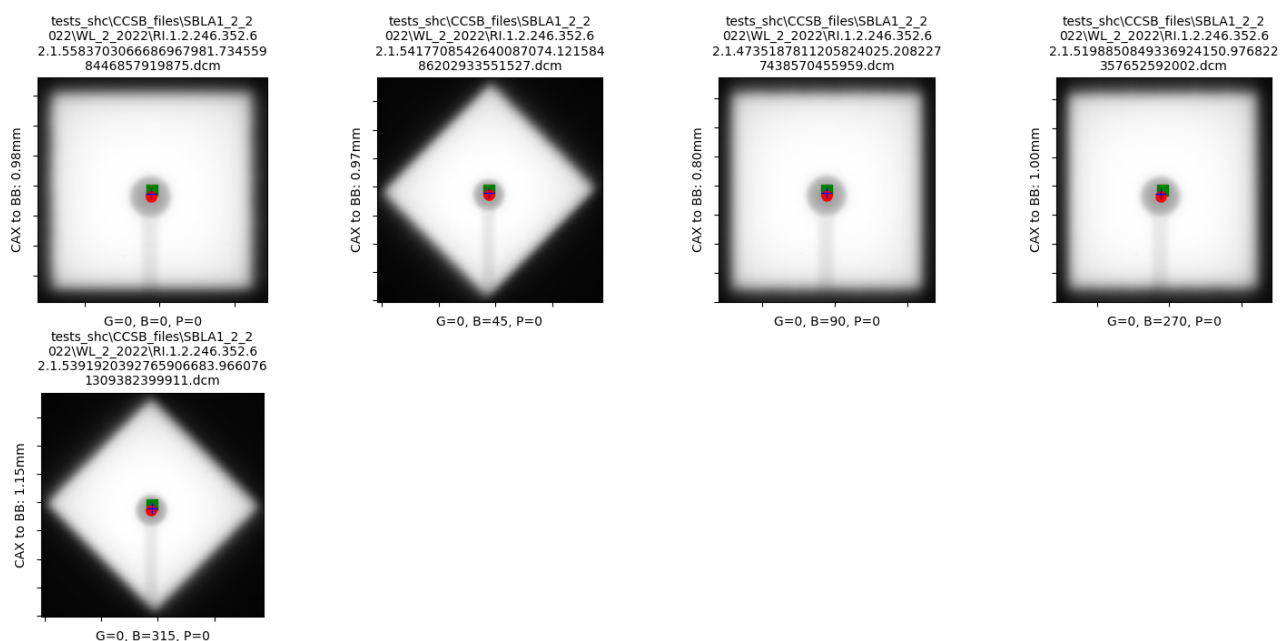

Couch images

