

Picket Fence results:

Magnification factor (SID/SAD): 1.50

Tolerance (mm): 0.5

Leaves passing (%): 100.0



Absolute median error (mm): 0.041

Mean picket spacing (mm): 14.9

Maximum error (mm): 0.219 on Picket 9, Leaf [10]

Gantry Angle: 187.04

Collimator Angle: 0.00

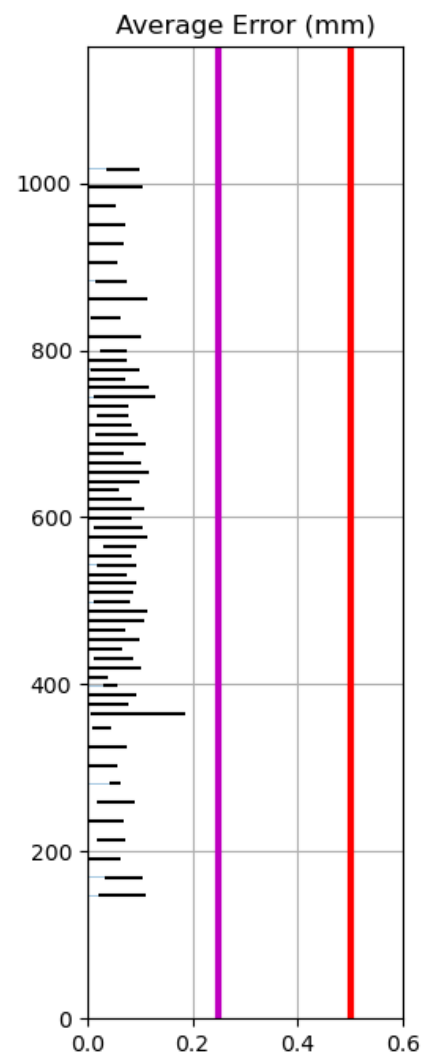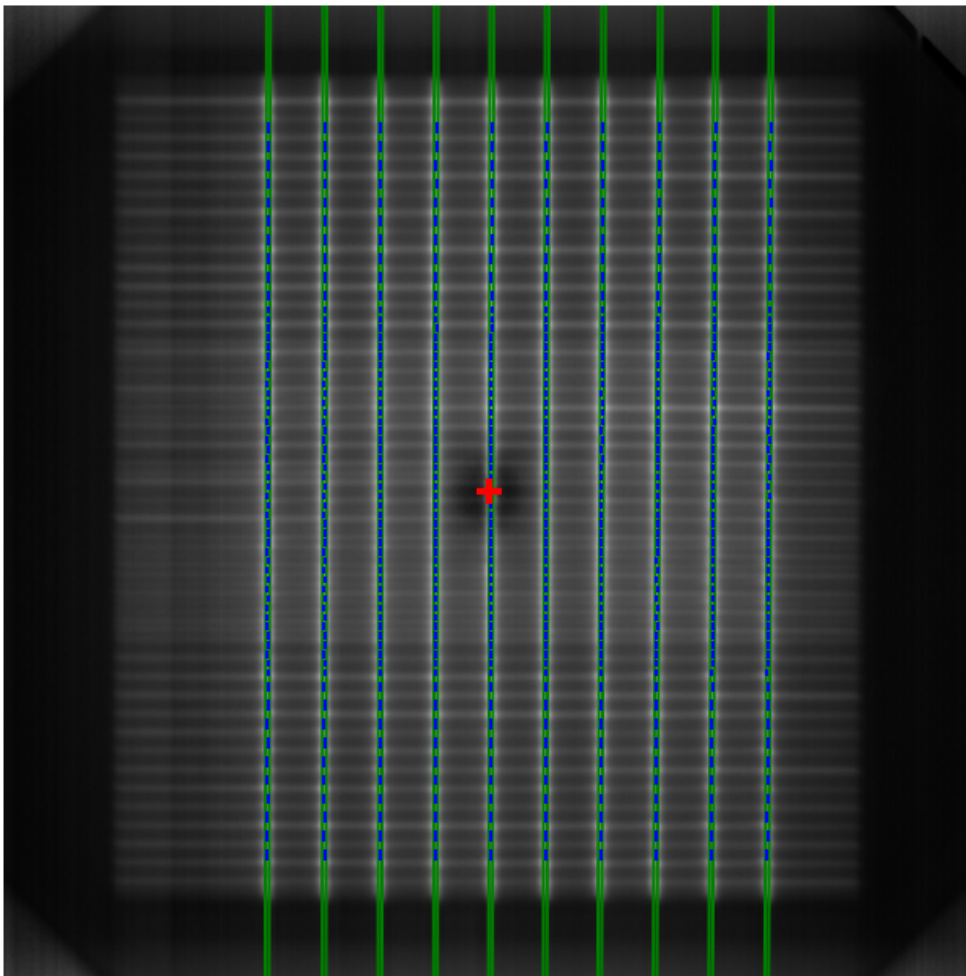

- CTP404 Results -



HU Linearity tolerance: 40

HU Linearity ROIs: Air: -1000.0, PMP: -193.0, LDPE: -97.0, Poly: -35.0, Acrylic: 132.0, Delrin: 376.0, Teflon: 1018.5

Geometric node spacing (mm): 50.03

Slice thickness (mm): 1.78

Low contrast visibility: 1.87

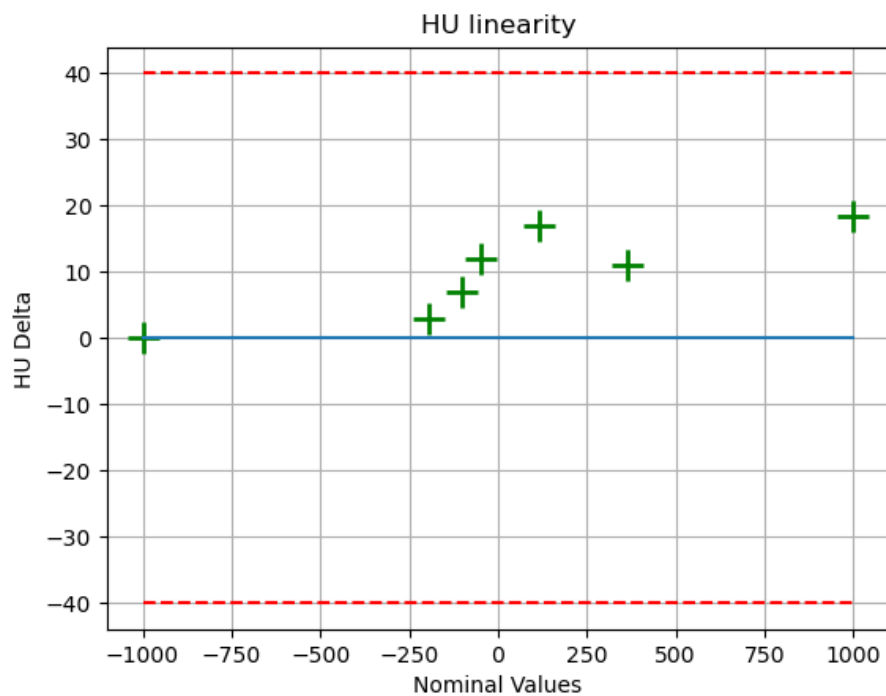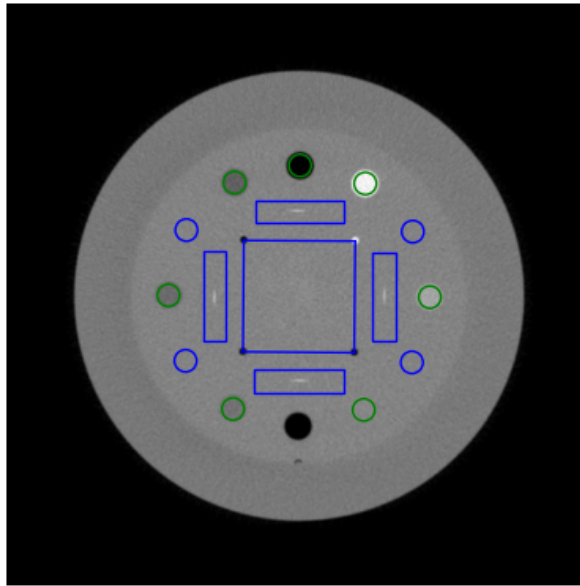



Notes:
testing

- CTP528 Results -

MTF 80% (lp/mm): 0.33

MTF 50% (lp/mm): 0.45

MTF 30% (lp/mm): 0.54

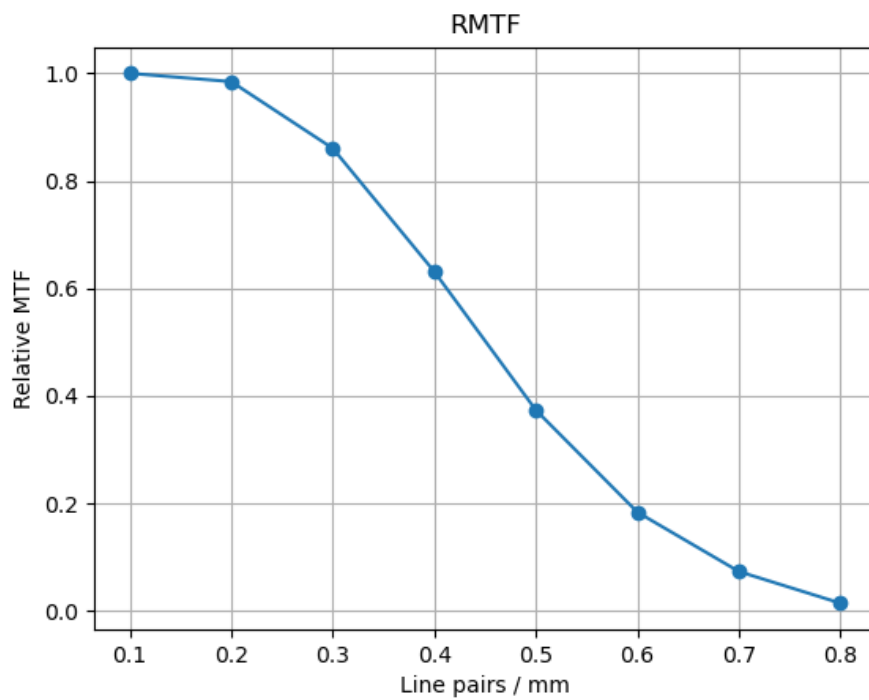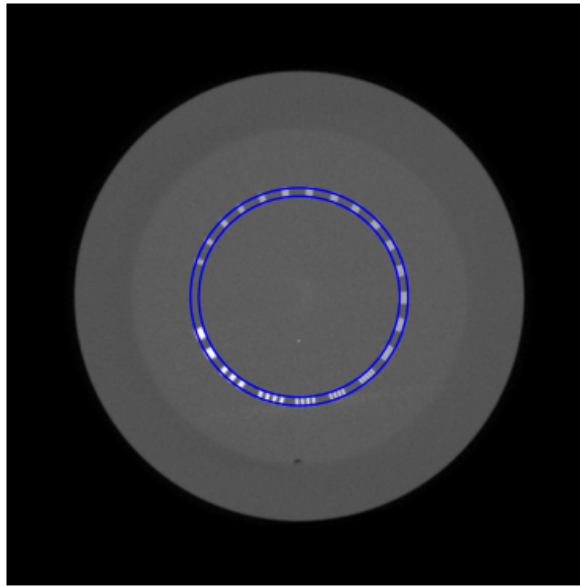

CatPhan 504 Analysis



- CTP486 Results -

Uniformity tolerance: 40


Uniformity ROIs: Top: 13.0, Right: 14.0, Bottom: 14.0, Left: 12.0, Center: 26.0

Uniformity Index: -1.36

Integral non-uniformity: 0.0069

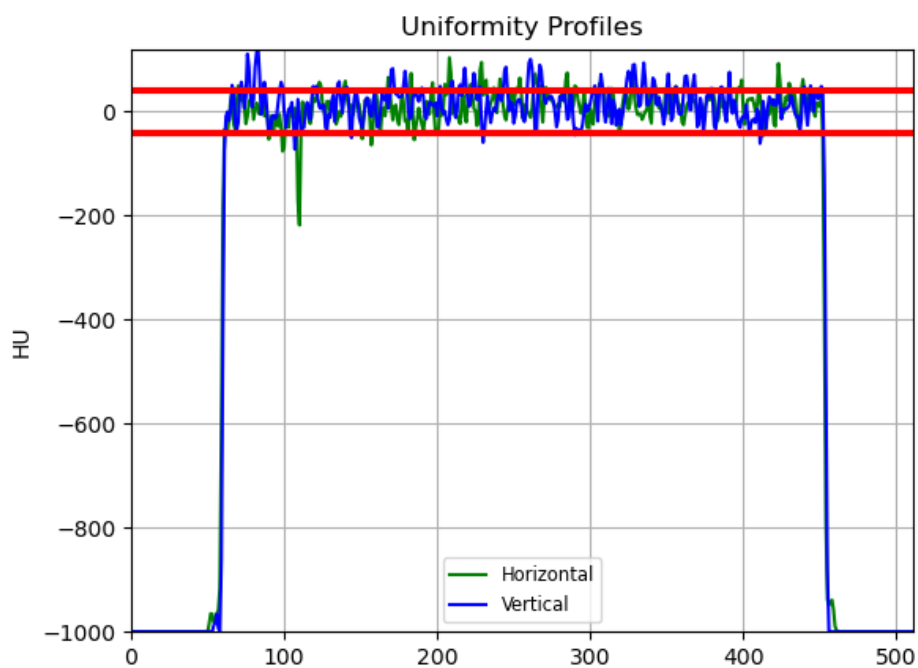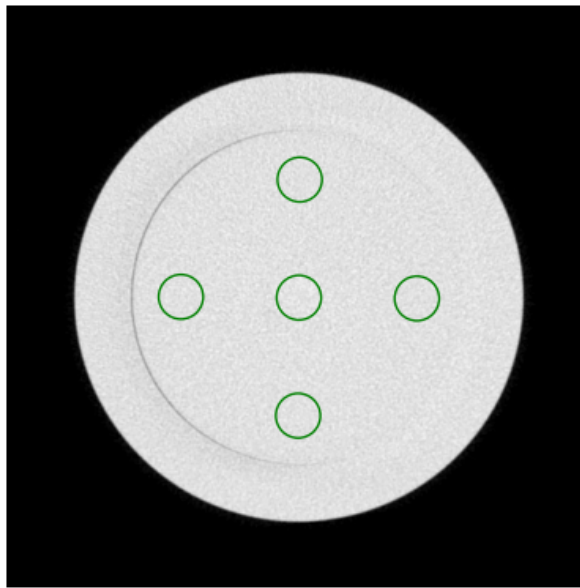

- CTP515 Results -
CNR threshold: 15
Low contrast ROIs "seen": 0

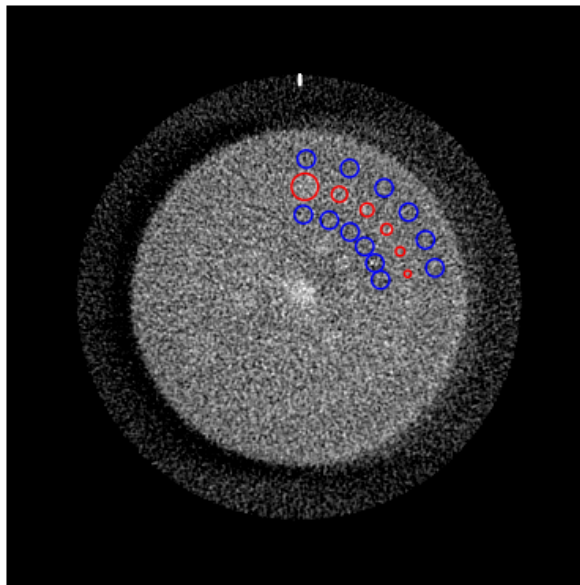

Dose Rate & Gantry Speed VMAT results:

Source-to-Image Distance (mm): 1500


Tolerance (%): 1.5

Absolute mean deviation (%): 0.62

Maximum deviation (%): 0.99


Open Image

20211223_181915_6x [MV]_G243_C360_T360_1.dcm

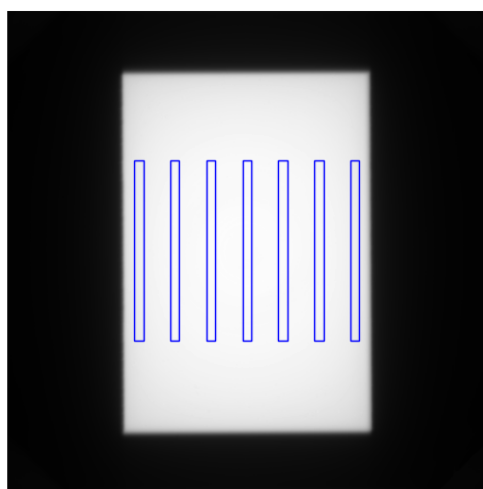


DMLC Image

20211223_181804_6x [MV]_G243_C360_T360_T3

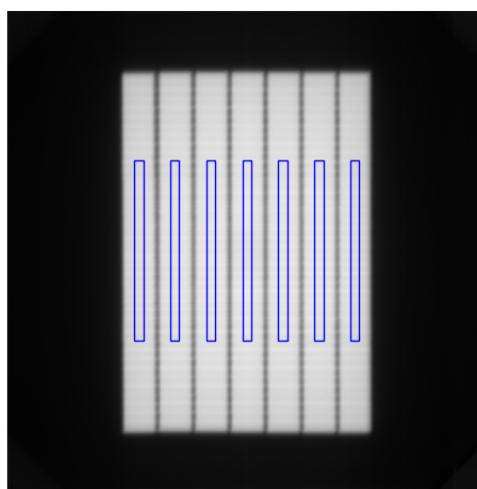

Median profiles

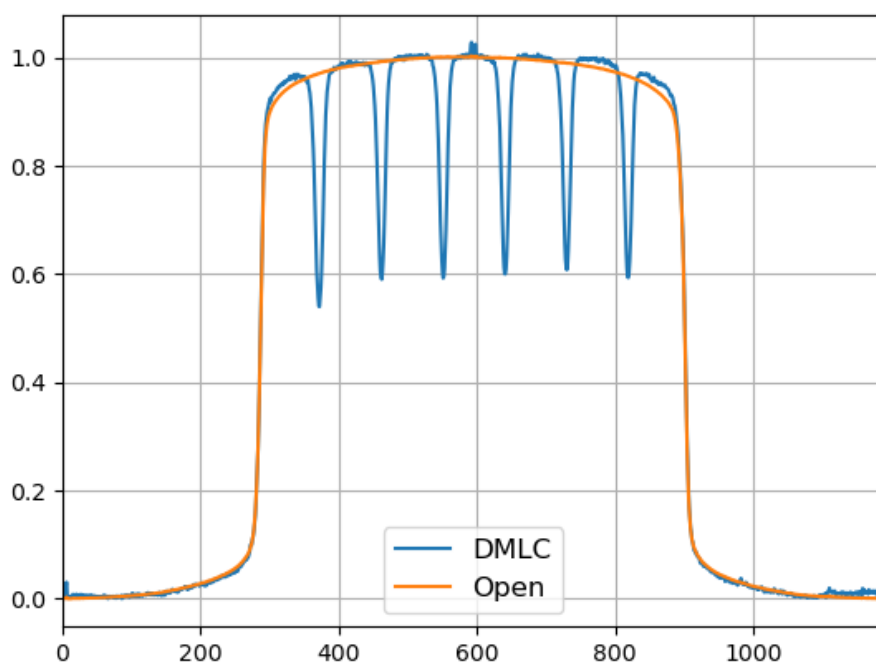

Dose Rate & MLC Speed VMAT results:

Source-to-Image Distance (mm): 1500


Tolerance (%): 1.5

Absolute mean deviation (%): 0.45

Maximum deviation (%): 0.89


Open Image

20220120_175430_6x [MV]_G32_C360_T0_2.dcm

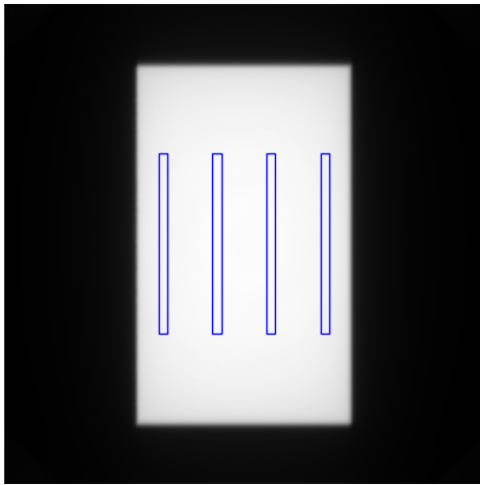


DMLC Image

20220120_175353_6x [MV]_G32_C360_T0_2.dcm

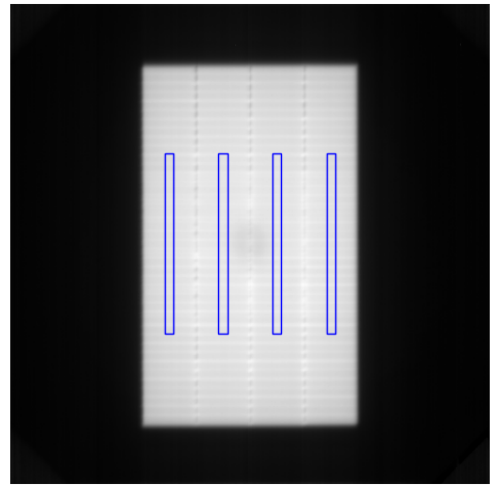

Median profiles







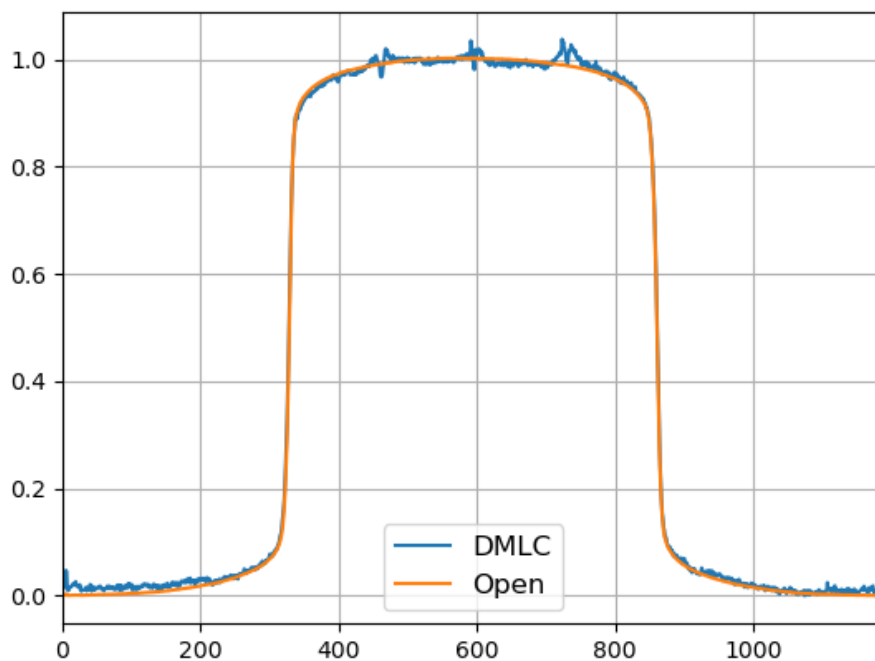


Winston-Lutz Analysis

Winston-Lutz Analysis


=====

Number of images: 12
Maximum 2D CAX->BB distance: 0.67mm
Median 2D CAX->BB distance: 0.31mm
Shift to iso: facing gantry, move BB: RIGHT 0.10mm; IN 0.25mm; UP 0.04mm
Gantry 3D isocenter diameter: 0.61mm (4/12 images considered)
Maximum Gantry RMS deviation (mm): 0.17mm
Maximum EPID RMS deviation (mm): 0.68mm
Gantry+Collimator 3D isocenter diameter: 0.69mm (8/12 images considered)
Collimator 2D isocenter diameter: 0.28mm (5/12 images considered)
Maximum Collimator RMS deviation (mm): 0.52
Couch 2D isocenter diameter: 0.62mm (5/12 images considered)
Maximum Couch RMS deviation (mm): 0.67

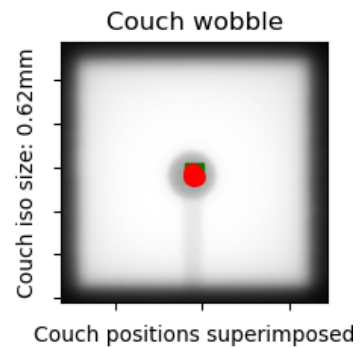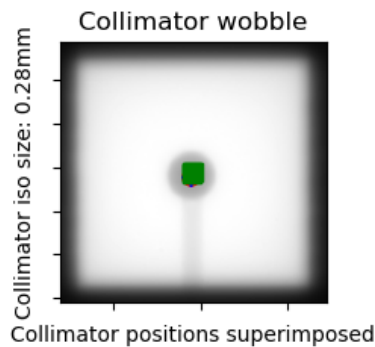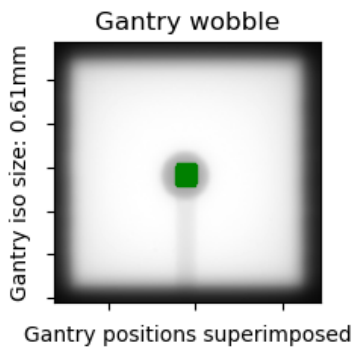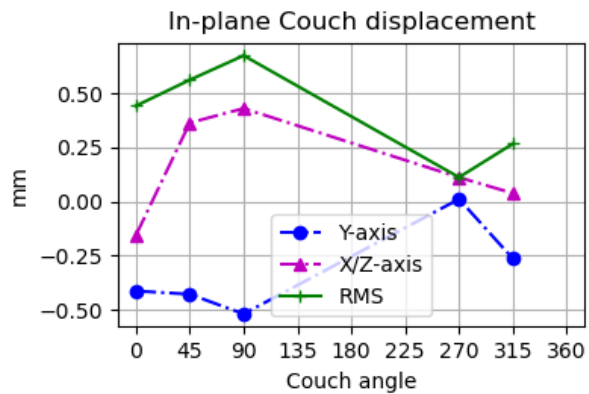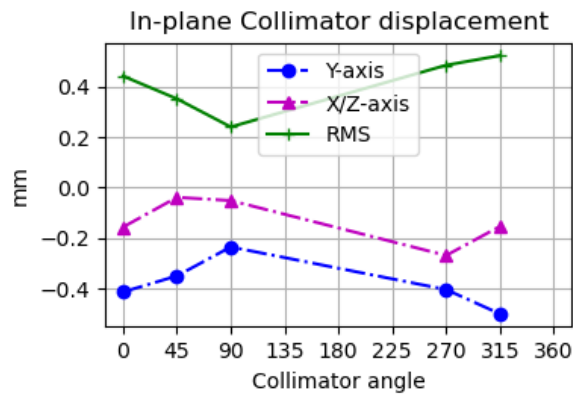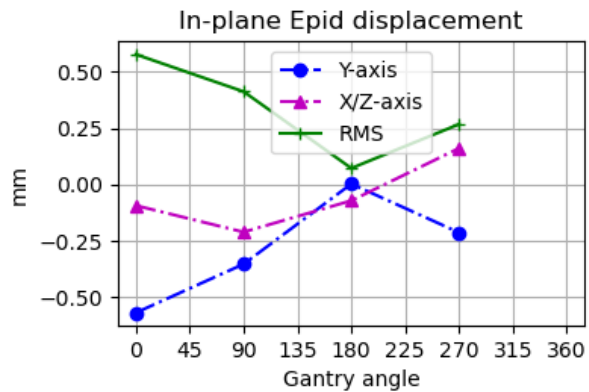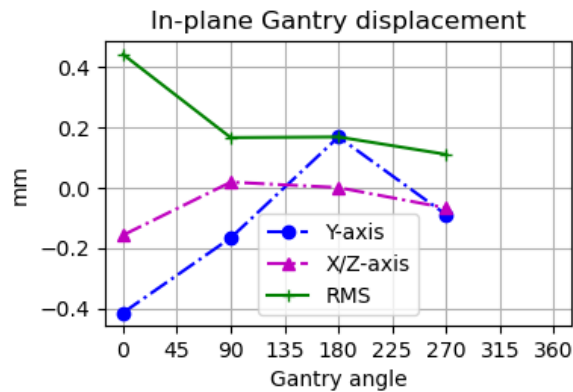


SI QC-3 results: ■ File: C:\Users\kjaps\Documents\Coding\GitHub\p...9_175938_2.5x [MV]


SI QC-3 Phantom Analysis

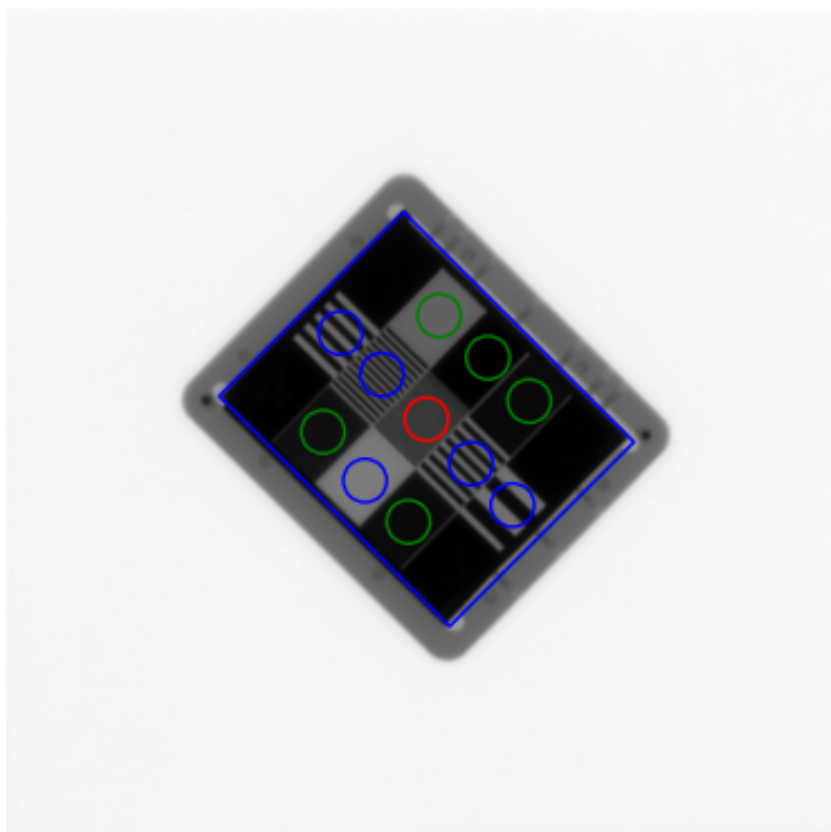



SI QC-kV results: ■ File: C:\Users\kjaps\Documents\Coding\GitHub\p...19_180233_100 [kV


SI QC-kV Phantom Analysis

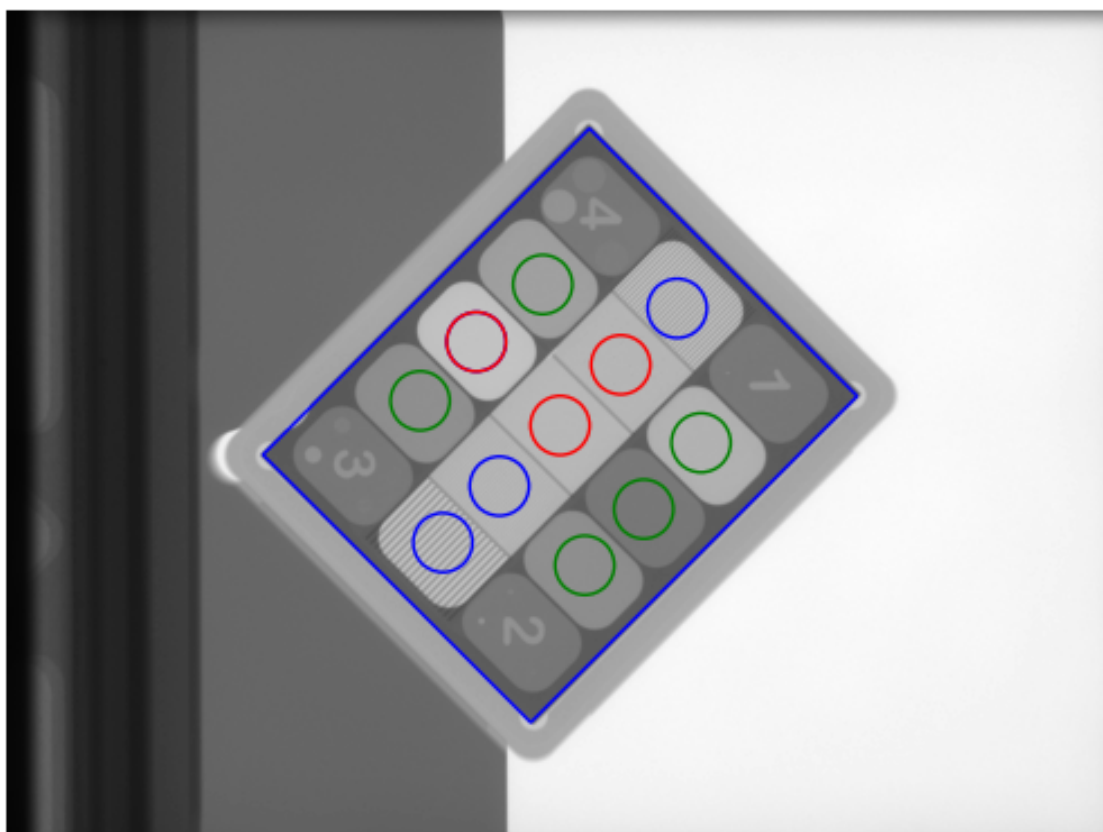

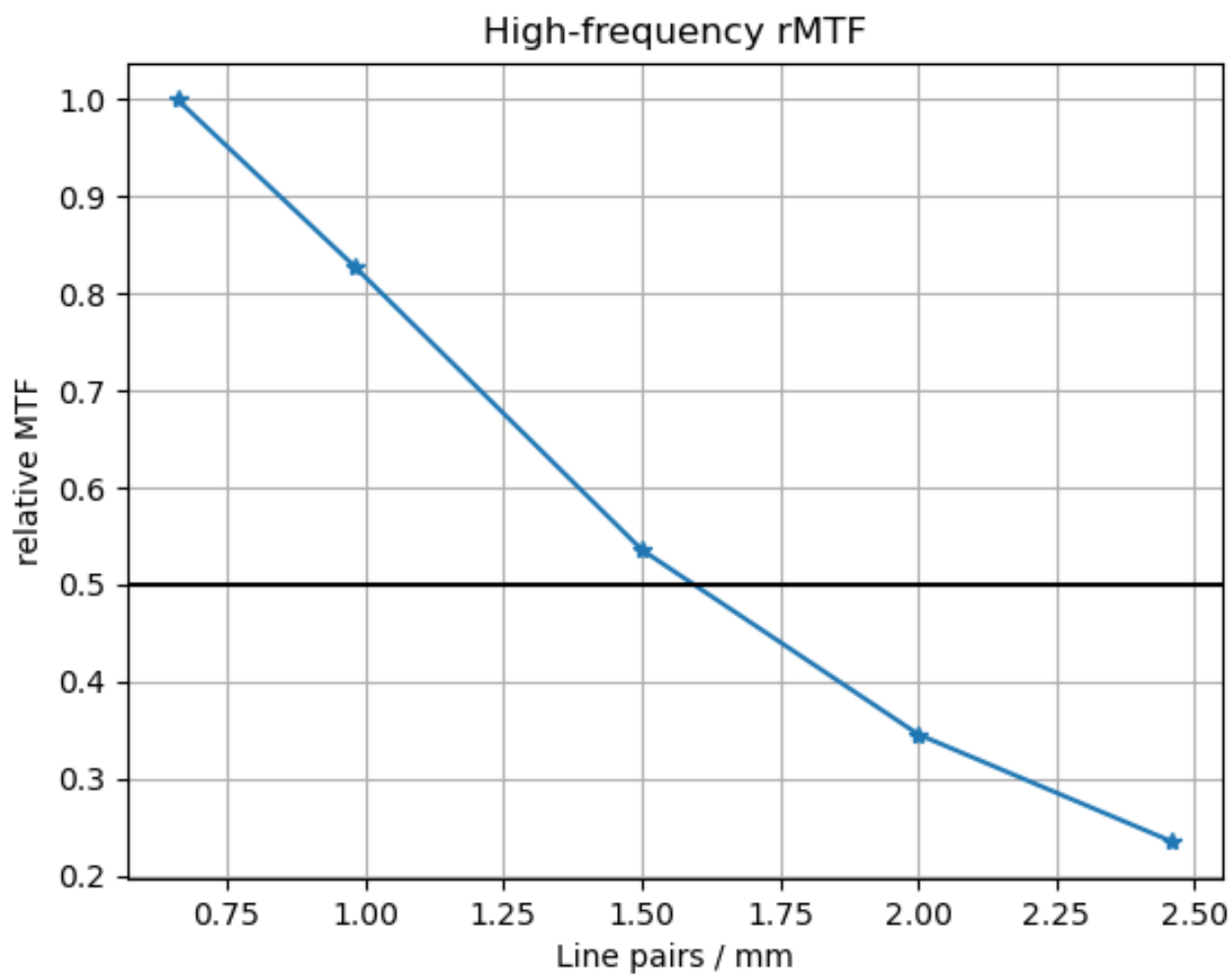


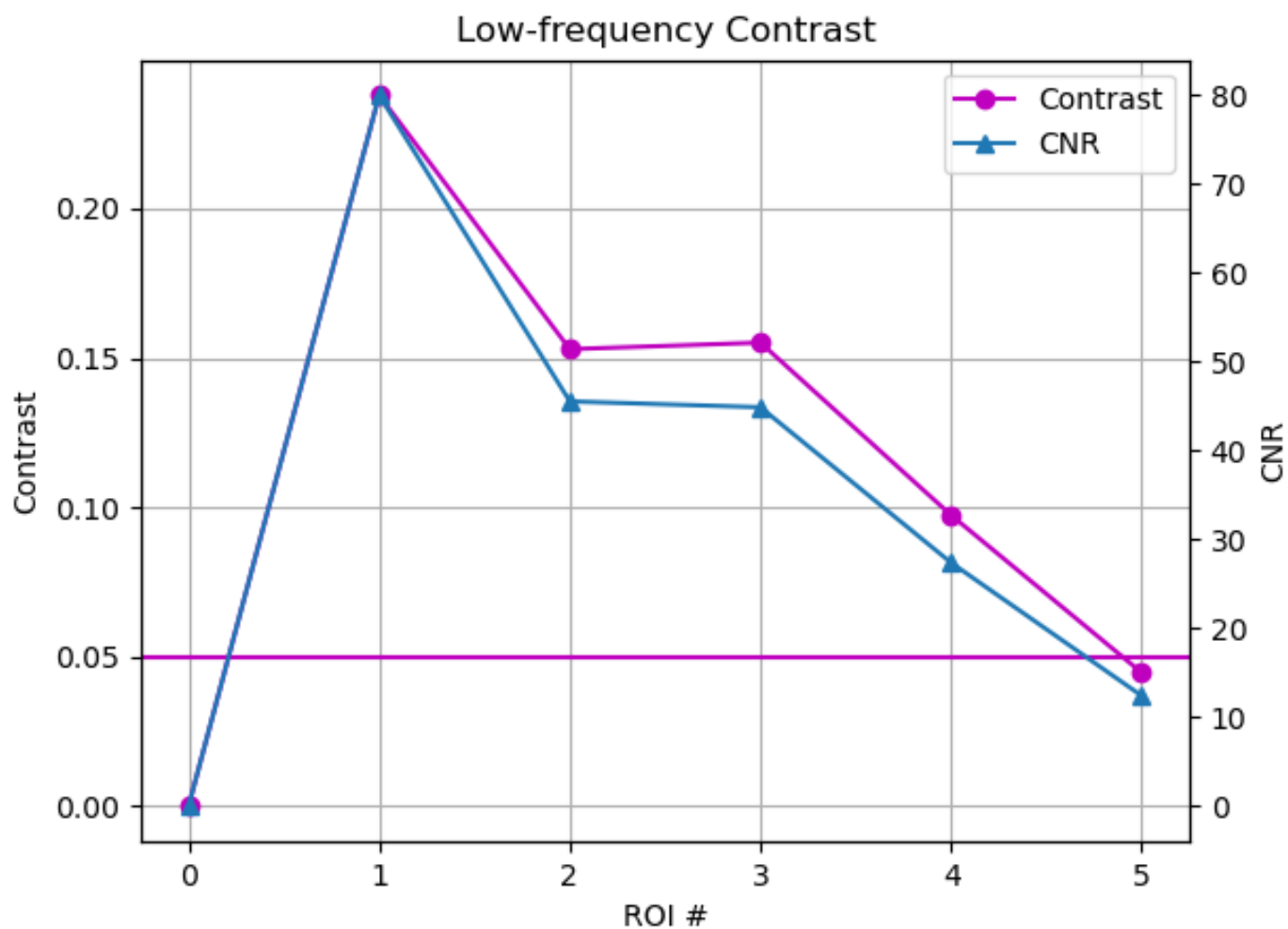

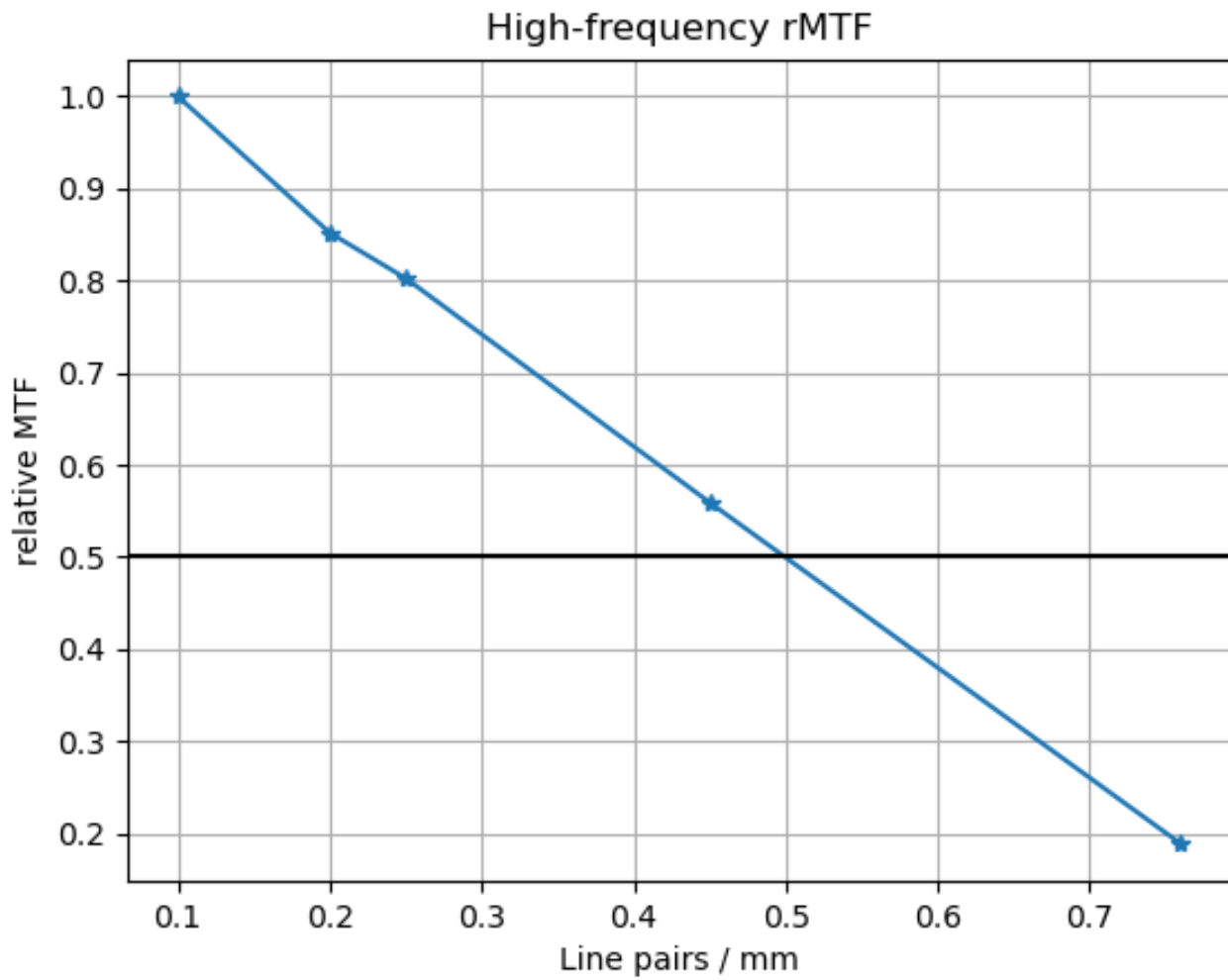


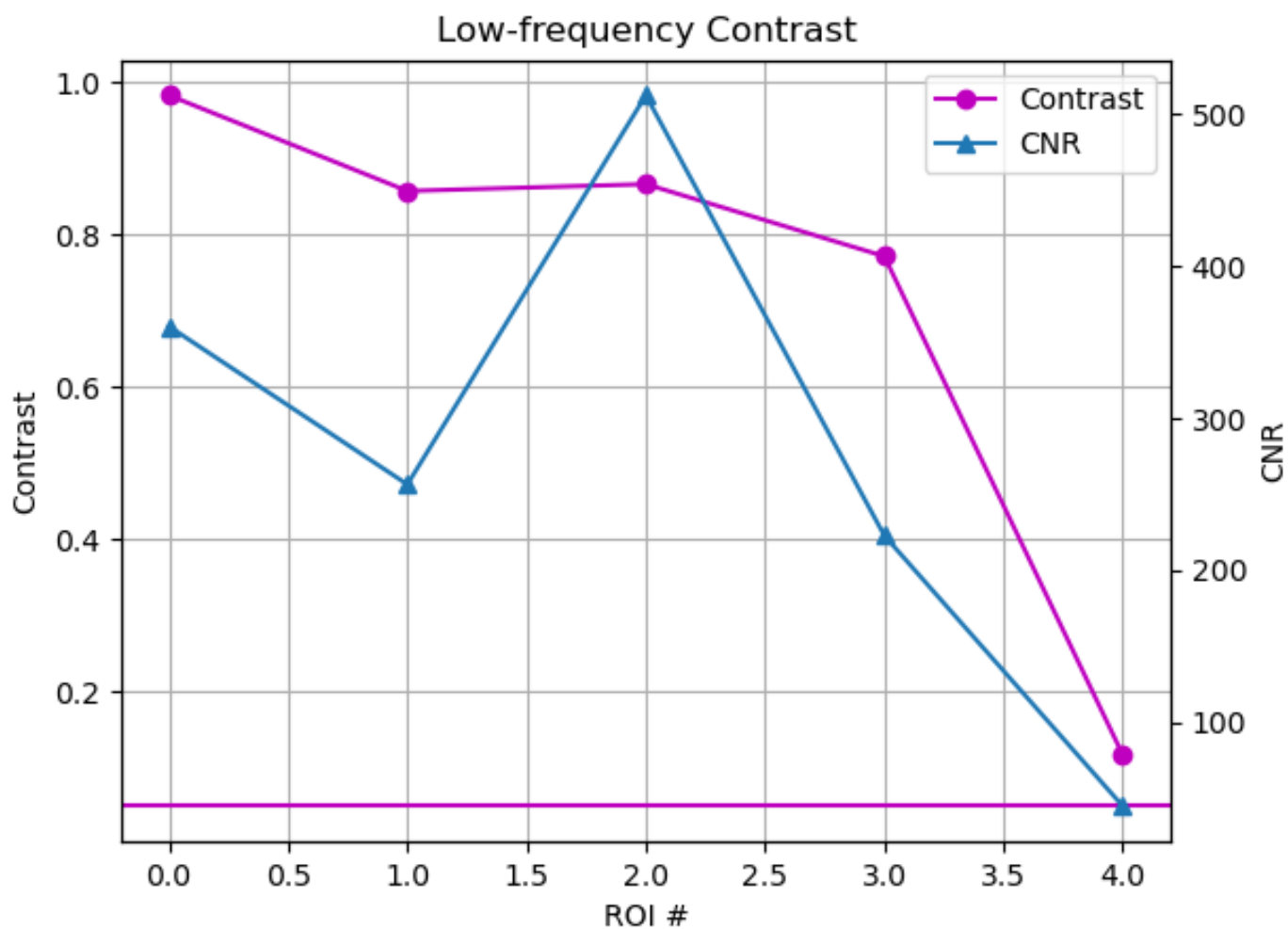

Gantry images

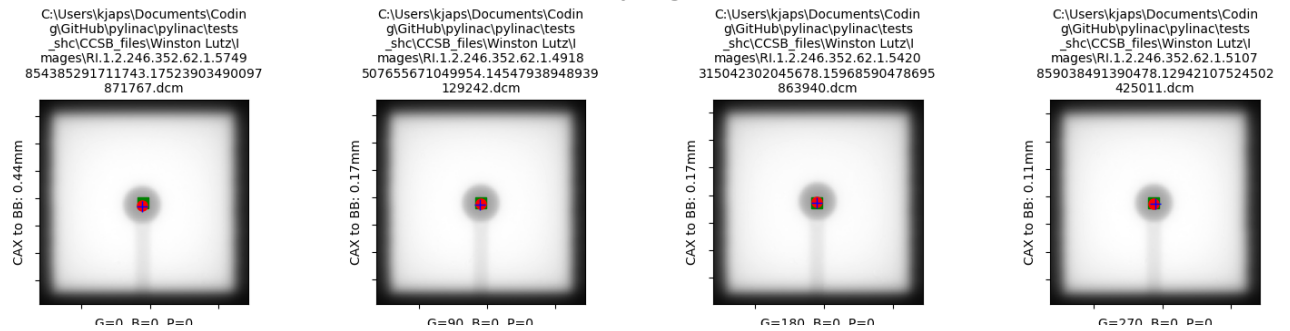

Collimator images

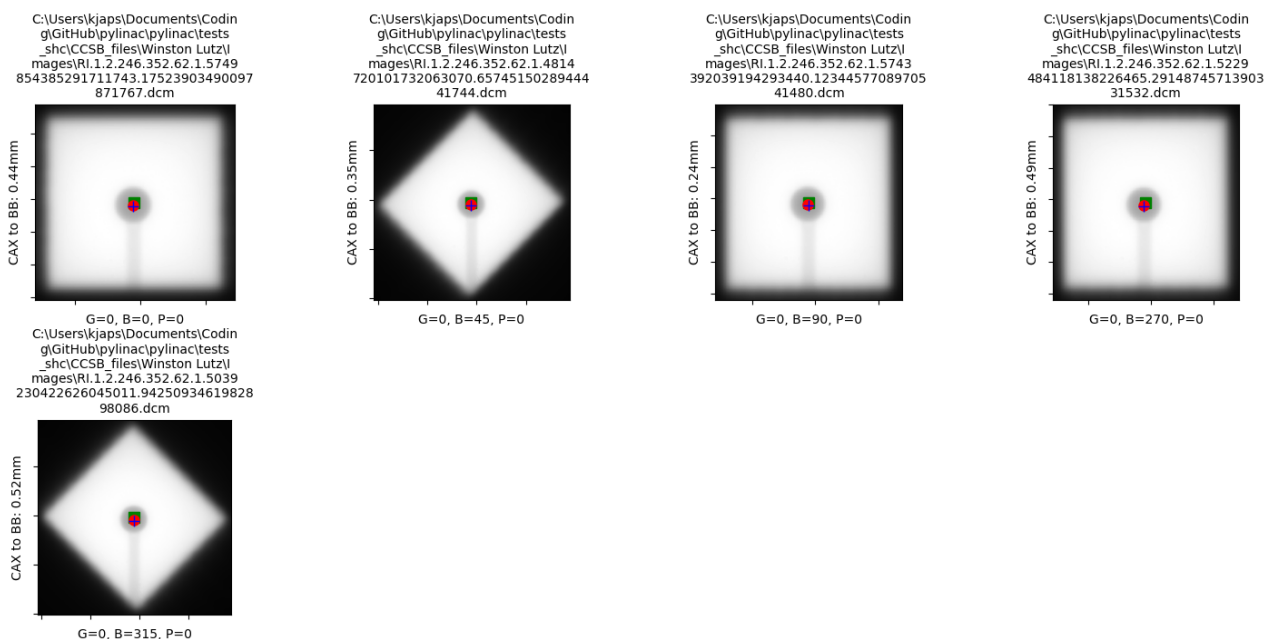

Couch images

